




LTP Director

New software for LTP and LTD studies

- Study long-term potentiation and long-term depression
- Control stimulation
- Control drug delivery and perfusion
- Acquire data
- Analyze online and offline

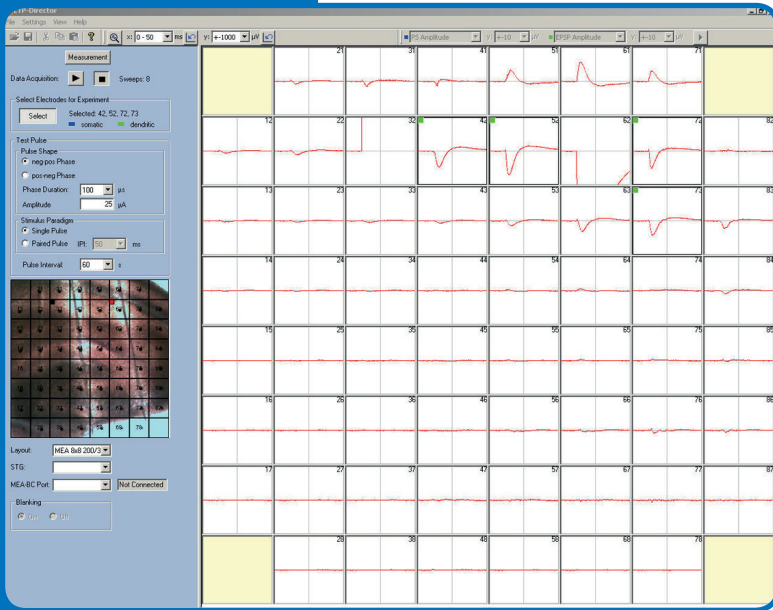

Test mode


In the test mode you choose all settings before you start recording:

- Load the image of your slice into the background
- Select stimulation electrodes for two independent pathways
- Select four electrodes for online analysis

Measurement mode

In the measurement mode you run the experiment and define all analysis parameters:

- Define your stimulation parameters
- Fill in the lab book
- Set up the experiment
- Define drug application
- Run experiment automatically
- Create your report sheet

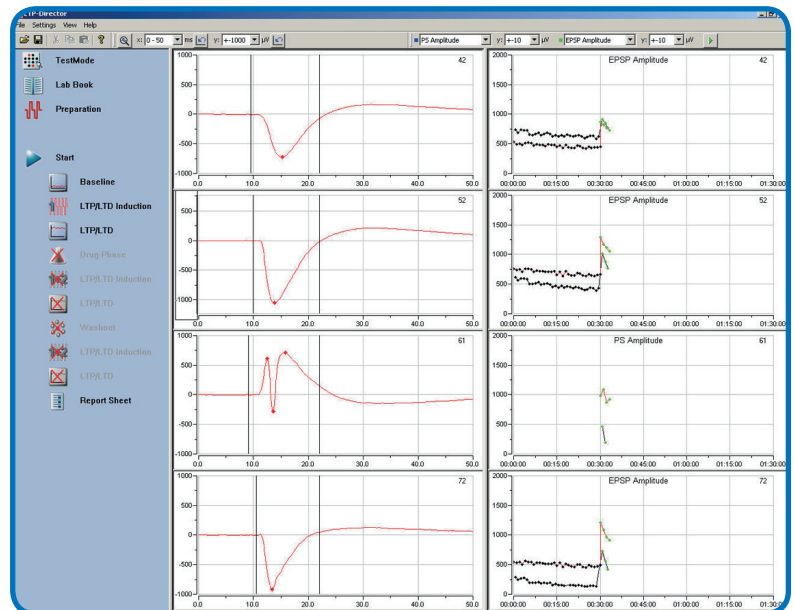

- Made
- in
- Germany

© 2013
Multi Channel Systems MCS GmbH

Product information is subject to change without notice. Products that are referred to in this document may be either trademarks and/or registered trademarks of the respective owners. The publisher and the author make no claim to these trademarks.