



## 9well-CC


### Culture Chamber for 256-9wellMEAs

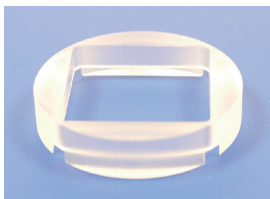

The 9well-CC culture chamber for 256-9wellMEAs is made of macrolon.

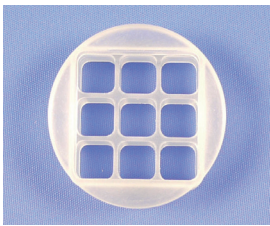

The 9well-CC culture chamber, mounted on a 9well macrolon quadrat.

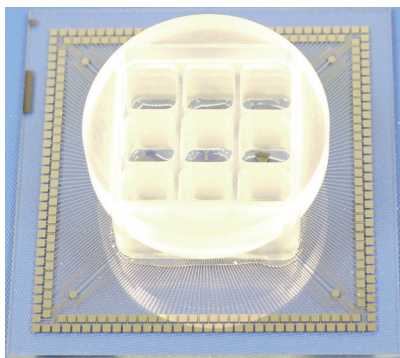

The 9well-CC mounted together with a 9well macrolon quadrat on a 256-9wellMEA300/30iR-ITO.

### Application

The 9well-CC culture chamber for 256-9wellMEA300/30iR-ITO is suitable for 9well macrolon quadrats. To use the macrolon quadrat in combination with the 9well-CC as a culture chamber, please insert a foil between macrolon quadrat and 9well-CC ring.

**Important:** The foil to cover the 9well macrolon quadrat with the 9well-CC has to be ordered separately!

For example: ALA MEA-MEM-SHEET (Set of 10 3" x 3" sheets of fluorinated ethylene-propylene membranes for sealed MEA culture dishes). Can be ordered at Multi Channel Systems MCS GmbH.

Scope of delivery: 1 9well-CC