


256 t MEA

Microelectrode Array with transparent Electrodes
for Use with MEA2100-256-System

256tMEA400/30iR-ITO

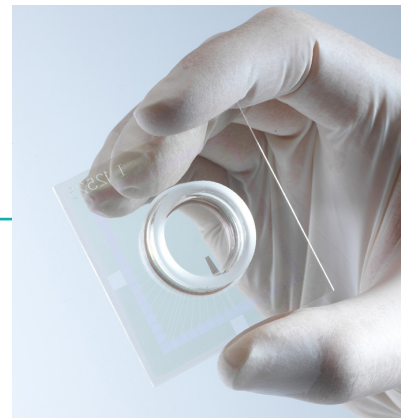


Technical Specification

General Characteristics	
Temperature compatibility	0 - 121 °C
Dimensions (W x D x H)	49 mm x 49 mm x 1 mm
60 t MEA Layout	
Base material	Glass
Electrode material	transparent TiN (Titanium nitride)
Electrode height	planar
Track material and contact pads	ITO (Indium tin oxide)
Electrode diameter	30 µm
Interelectrode distance (center to center)	400 µm
Isolation material	Silicon nitride (SiN)
Electrode impedance	< 250 kΩ
Electrode layout grid	16 x 16
Number of recording electrodes	252
Number of reference electrodes	4 internal reference electrode (iR)
Contact pads for reference electrodes, connected to ground	4
Software: MC_Rack	MC_Rack: 2 dim. (MEA) or Configuration
Software: Multi Channel Suite	Multi Channel Experimenter: MEA Configuration
MEA Perfusion Chamber	
(w/o) Without ring (gr) Glass ring Note: On explicit customers order only! The glue for the glass ring may cover the reference electrodes, which make them inefficient. You have to use an AgAgCl pellet as reference. (pr) Plastic ring without thread (pr-T) Plastic ring with thread	ID +/- 19 mm, OD +/- 24 mm, height 6 / 12 mm ID +/- 26.5 mm, OD +/- 30 mm, height 6 / 12 mm ID +/- 26 mm, OD +/- 30 mm, height 6 / 12 mm
Sterilization	These MEA types are heatstable up to 121 °C! They can be autoclaved and coated with different procedures for cell and tissue cultures.

256 t MEA400/30iR-ITO: Electrode Layout

Electrode Layout

Important: Please do not treat a MEA with transparent electrodes in a Plasma Cleaner!
 To make it hydrophilic use PBS over night.

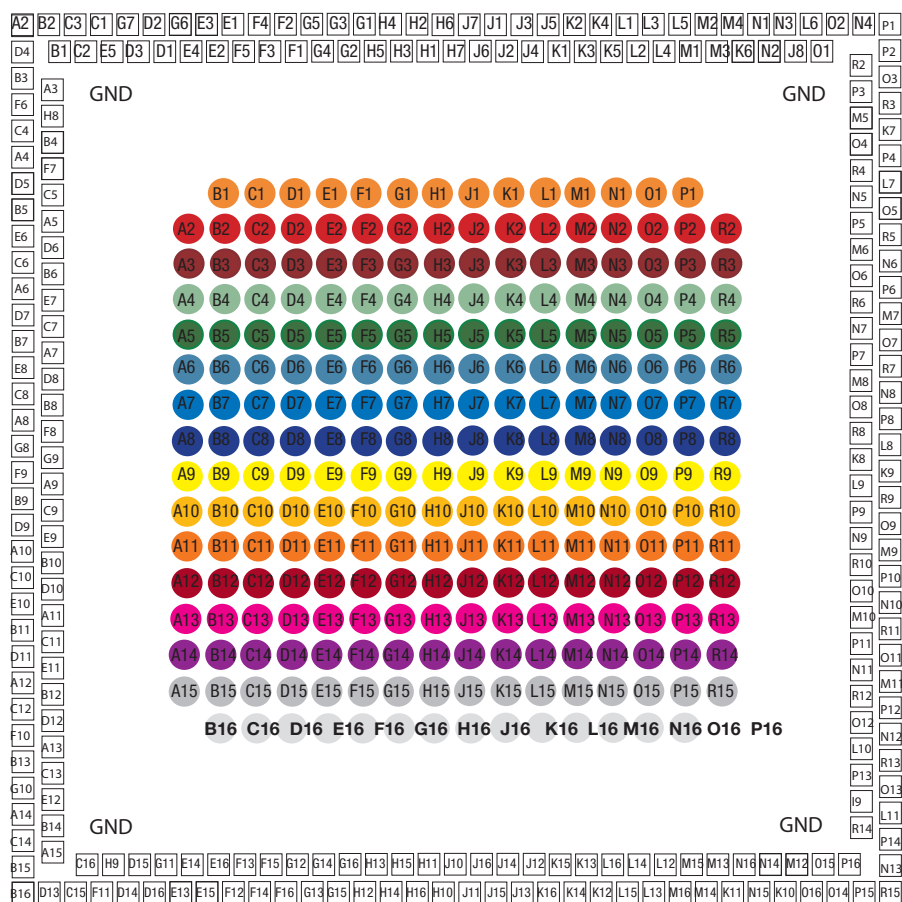

Technical Specifications

Numbering	The numbering of 256tMEA electrodes in the 16 x 16 grid follows the standard numbering scheme for square grids: The number code is the hardware identifier of the electrode. The reference electrodes (REF) are connected to ground (GND).
256MEAs are symmetrical	The 256MEA is the only MEA type, which is rotationally symmetrical. So, the orientation does not matter.

256 t MEA: Table of Hardware ID and Electrode

Connector Electrode ID	Socket 1 Stimulation Socket	Spring Contact	Connector Electrode ID	Socket 2 Stimulation Socket	Spring Contact	Connector Electrode ID	Socket 3 Stimulation Socket	Spring Contact	Connector Electrode ID	Socket 4 Stimulation Socket	Spring Contact
A2	2	1	J 9	60	191	C15	61	98	A10	25	113
B1	4	133	K7	7	38	C16	62	222	A11	22	232
B2	1	2	K8	34	178	D13	63	99	A12	15	108
C1	5	4	K9	33	51	D14	57	96	A13	12	227
C2	6	134	L10	56	189	D15	58	220	A14	5	103
C3	3	3	L11	59	64	D16	55	95	A15	4	223
D1	12	137	L7	11	40	E13	53	94	A3	62	252
D2	9	6	L8	31	50	E14	54	218	A4	55	128
D3	10	136	L9	36	179	E15	51	93	A5	52	247
E1	15	9	M10	46	184	E16	52	217	A6	45	123
E2	16	139	M11	49	59	F11	59	97	A7	42	242
E3	13	8	M5	8	165	F12	49	92	A8	35	118
E4	14	138	M6	18	170	F13	50	216	A9	32	237
E5	8	135	M7	21	45	F14	47	91	B10	26	234
F1	22	142	M8	28	175	F15	48	215	B11	19	110
F2	19	11	M9	39	54	F16	45	90	B12	16	229
F3	20	141	N10	43	56	G11	56	219	B13	9	105
F4	17	10	N11	50	186	G12	46	214	B14	6	224
F5	18	140	N12	53	61	G13	43	89	B15	1	101
G1	25	14	N13	63	66	G14	44	213	B16	2	100
G2	26	144	N5	14	168	G15	41	88	B3	61	131
G3	23	13	N6	17	43	G16	42	212	B4	58	250
G4	24	143	N7	24	173	H10	33	84	B5	51	126
G5	21	12	N8	27	48	H11	36	209	B6	48	245
G6	11	7	N9	40	181	H12	39	87	B7	41	121
G7	7	5	O10	44	183	H13	40	211	B8	38	240
H1	32	147	O11	47	58	H14	37	86	B9	29	115
H2	29	16	O12	54	188	H15	38	210	C10	23	112
H3	30	146	O13	57	63	H16	35	85	C11	20	231
H4	27	15	O3	3	36	H9	60	221	C12	13	107
H5	28	145	O4	10	166	J10	34	208	C13	10	226
H6	31	17	O5	13	41	J11	31	83	C14	3	102
H7	34	148	O6	20	171	J12	28	205	C4	57	129
J 1	35	19	O7	23	46	J13	27	81	C5	54	248
J 2	38	150	O8	30	176	J14	30	206	C6	47	124
J 3	37	20	O9	37	53	J15	29	82	C7	44	243
J 4	40	151	P1	2	34	J16	32	207	C8	37	119
J 5	39	21	P10	41	55	K10	7	71	C9	30	236
J 6	36	149	P11	48	185	K11	11	73	D10	24	233
J 7	33	18	P12	51	60	K12	21	78	D11	17	109
J 8	60	161	P13	58	190	K13	24	203	D12	14	228
K1	42	152	P14	61	65	K14	23	79	D4	63	132
K2	41	22	P2	1	35	K15	26	204	D5	53	127
K3	44	153	P3	6	164	K16	25	80	D6	50	246
K4	43	23	P4	9	39	L12	18	200	D7	43	122
K5	46	154	P5	16	169	L13	17	76	D8	40	241
K6	56	159	P6	19	44	L14	20	201	D9	27	114
L1	45	24	P7	26	174	L15	19	77	E10	21	111
L2	48	155	P8	29	49	L16	22	202	E11	18	230
L3	47	25	P9	38	180	M12	8	195	E12	8	225
L4	50	156	R10	42	182	M13	14	198	E6	49	125
L5	49	26	R11	45	57	M14	13	74	E7	46	244
L6	59	31	R12	52	187	M15	16	199	E8	39	120
M1	52	157	R13	55	62	M16	15	75	E9	28	235
M2	51	27	R14	62	192	N14	10	196	F10	11	106
M3	54	158	R2	4	163	N15	9	72	F6	59	130
M4	53	28	R3	5	37	N16	12	197	F7	56	249
N1	55	29	R4	12	167	O14	3	69	F8	36	239
N2	58	160	R5	15	42	O15	6	194	F9	31	116
N3	57	30	R6	22	172	O16	5	70	G10	7	104
N4	63	33	R7	25	47	P15	1	68	G8	33	117
O1	62	162	R8	32	177	P16	4	193	G9	34	238
O2	61	32	R9	35	52	R15	2	67	H8	60	251
GND	64	253	GND	64	254	GND	64	255	GND	64	256



Multi Channel Systems
Multi Channel Systems
MCS GmbH
Aspenhaustrasse 21
72770 Reutlingen Germany

Sales:
sales@multichannelsystems.com
Technical Support:
support@multichannelsystems.com
Web:
multichannelsystems.com

Americas
Tel: (+1) 833 668 8632
Europe, Middle East, Africa
Tel: (+49) 7121 909 2525
Asia Pacific
Tel: (+86) 21 6226 0239

Copyright © 2025 Multi Channel Systems

Product information is subject to change without notice. Multi Channel Systems is a trademark of Harvard Bioscience, Inc. or its affiliated companies. Harvard is a registered trademark of Harvard University. The mark Harvard Bioscience is being used pursuant to a license agreement between Harvard University and Harvard Bioscience, Inc.

February 2025